


Oxigel Powder


OPEN SURGERY AND LAPAROSCOPIC APPLICATORS

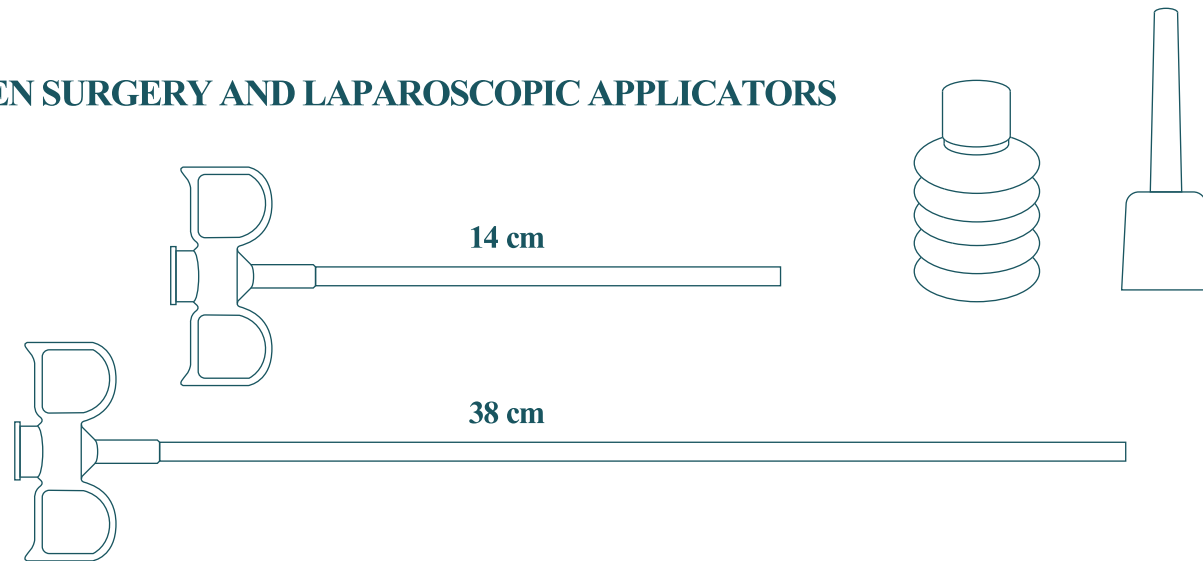


Λυκούργου 2, Νέα Ιωνία
Τηλ.: 216 8090015-16-17
Fax: 216 8090018
www.pelmyplusmedical.com
info@pelmyplusmedical.com

Betatech Medikal Cihazlar
San. Müm. İç ve Dış
Tic. Ltd.Şti.

BETATECH
MEDICAL.COM

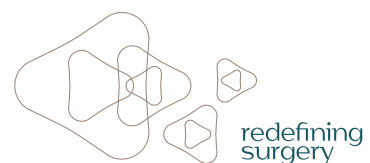

Oxigel Powder

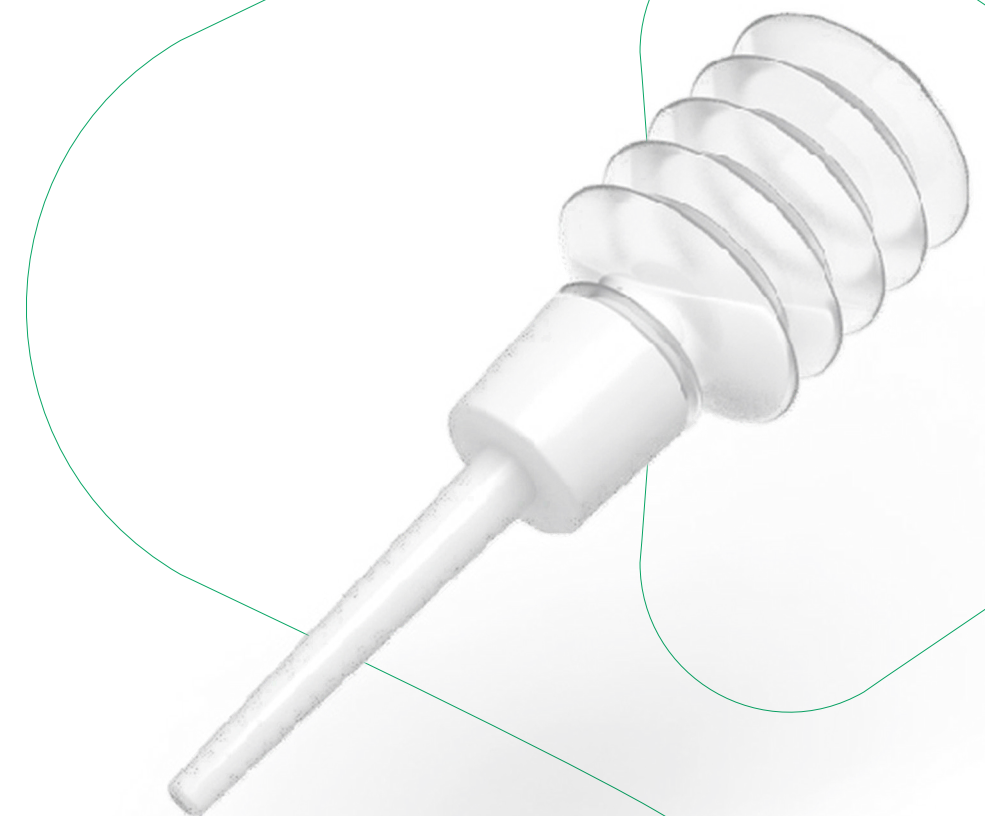

Latest Generation
Powder Hemostat



Oxcel Powder

The new generation of topical absorbable hemostats, Oxcel Powder, has been developed to control persistent oozing bleeding on large surfaces. The particles of Oxcel Powder are biocompatible, non-pyrogenic and derived from purified plant starch, so there are no human or animal materials.

Application Technique

1 DRY

Gauze dry the wound and remove all excess blood.

2 APPLY

Immediately after removing excess blood, apply OXICEL POWDER over the entire wound.

3 COMPRESS

Apply direct pressure over OXICEL POWDER using moist gauze for 1 to 2 minutes.

4 IRRIGATE

Remove excess OXICEL POWDER after hemostasis is achieved.

Immediate Hemostasis Absorption within hours

Oxcel Powder upon contact with blood, initiates the clotting process. Blood solids such as platelets, red blood cells and blood proteins are concentrated in powder particles to create a gelled matrix. The normal clotting mechanism is improved, regardless of the coagulation status of the patient, by creating a barrier to further blood loss. Oxcel Powder absorbed within 24-48 hours.

Proven protection and effectiveness in a range of surgical areas and types of procedures

- General surgery
- Gynecological
- Urology
- Cardiac and vascular surgery
- Neurosurgery
- Traumatology
- Orthopedicsurgery
- Plasticsurgery
- ENTsurgery



Safe

- 100% plant-based polysaccharides
- Naturally Safe
- Biocompatible
- Absorbable in 24-48 hours

Simple

- No additional prep time necessary
- No special storage conditions required
- No mixing or refrigeration
- Easy to use applicators

Effective

- Complete hemostasis in minutes
- Ensure wide area coverage on rough surfaces and in areas hard to reach
- Cost Effective

Delivery System:

Applying Oxcel Powder correctly and directly is quick and easy with the spray applicator. A light, plastic device with a long, flexible tube is featured in the distribution system. Extended reach applicator tips from Oxcel Powder Applicators provide accurate and precise delivery of Oxcel Powder hemostatic powder from a simple, single-use device.